


55" 3D 4K Medical Display

FM-A5505DGC

Rev. 02



- *OLED technology. High contrast ratio. Perfect blacks.*
- *Supports HDR for an incredible image.*
- *Immersive 3D viewing is ideally suited for large screens.*

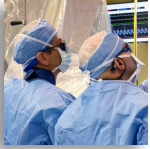

More Surface Textures More Depth of Field

FM-A5505DGC Rev.02, the 55" 3D display monitor from FSN, delivers OLED and HDR technology for true blacks, better contrast, power efficiency, and lighter weight. The immersive viewing experience of 3D is ideally suited for a large screen of this size.

FM-A5505DGC Rev.02 is compatible with highly specialized visualization equipment used in operating rooms and surgical suites. It is suitable for showing images at sizes big enough to remain visible across the operating room. Picture quality is crystal clear. The HDR wide color gamut is typically larger than standard color values used in Rec.709. DCI-P3 color space can reach 98%.




3D eye wear included.




Complete Control

A detailed user menu offers control for the best possible image settings. The interface will adjust Picture, Color, Advanced, Setup, Preset, Layout, and 3D.

FSN's attention to engineering details, such as rounded bezel corners, backlit touch buttons, and a comprehensive set of included accessories, allows for safe and efficient interaction. This large screen delivers everything that professionals need in a 3D, 4K UHD medical-grade monitor.

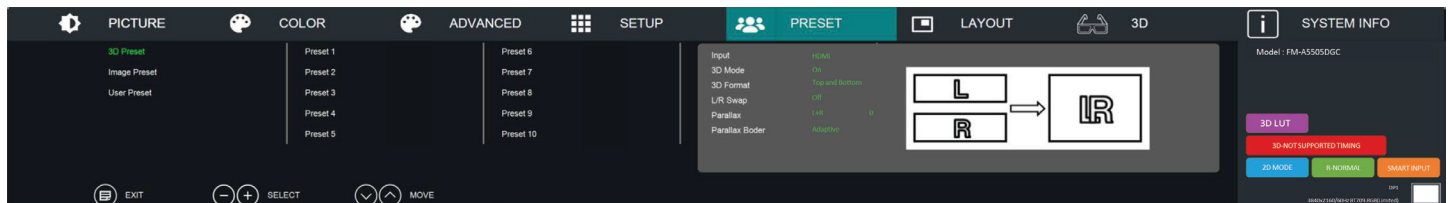

The FM-A5505DGC Rev.02 on-screen display (OSD) offers detailed user control menus.

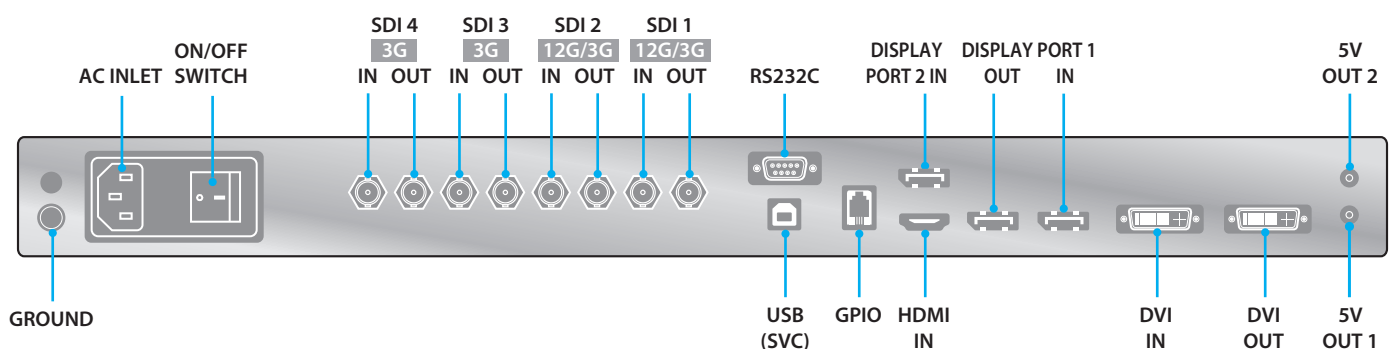

55" 3D 4K Medical Display

FM-A5505DGC

Rev. 02




Specifications

General Features

Item	Description												
Panel	55 inch OLED												
Resolution	3840 x 2160 pixels												
Aspect Ratio	16 : 9												
Active Area	1209.60 (H)mm x 680.40 (V)mm												
Pixel Pitch (mm)	0.315 x 0.315												
Response Time (typical)	1 ms (gray to gray)												
Number of Colors	1.07 Billion												
Luminance, Bypass (2D / 3D) cd/m ²	<table><tr><td></td><td><u>Minimum</u></td><td><u>Typical</u></td><td><u>Measurement</u></td></tr><tr><td>Peak</td><td>376/151</td><td>470/192</td><td>25% white box</td></tr><tr><td>Normal</td><td>144/58</td><td>80/73</td><td>100% white pattern</td></tr></table>		<u>Minimum</u>	<u>Typical</u>	<u>Measurement</u>	Peak	376/151	470/192	25% white box	Normal	144/58	80/73	100% white pattern
	<u>Minimum</u>	<u>Typical</u>	<u>Measurement</u>										
Peak	376/151	470/192	25% white box										
Normal	144/58	80/73	100% white pattern										
Luminance, D65 (2D / 3D) cd/m ²	<table><tr><td>Peak</td><td>366/135</td><td>420/174</td><td>25% white box</td></tr><tr><td>Normal</td><td>136/55</td><td>170/69</td><td>100% white pattern</td></tr></table>	Peak	366/135	420/174	25% white box	Normal	136/55	170/69	100% white pattern				
Peak	366/135	420/174	25% white box										
Normal	136/55	170/69	100% white pattern										
Luminance, D65 (2D / 3D, HDR) cd/m ²	<table><tr><td>Peak</td><td>496/202</td><td>620/252</td><td>10% white box</td></tr><tr><td>Normal</td><td>336/135</td><td>420/170</td><td>25% white pattern</td></tr></table>	Peak	496/202	620/252	10% white box	Normal	336/135	420/170	25% white pattern				
Peak	496/202	620/252	10% white box										
Normal	336/135	420/170	25% white pattern										
Gamut	BT.709 and BT.2020 compatible												
Contrast Ratio, D65 (2D, typical) cd/m ²	420,000 : 1 (peak), 170,000 : 1 (normal)												
Surface Treatment	Anti-glare												
Viewing Angle	(2D) R/L 120°, U/D 120° (3D) U/D 16°												

FSN Medical Technologies products meet EU RoHS Directive 2011/65/EU+2015/863/EU and Registration, Evaluation, Authorisation and Restriction of Chemicals (REACH) obligations. FSN manages the information flow as required by regulations.

Item	Description
Input Signal	1 x HDMI (2.0, HDCP 2.3) 2 x DP (1.2 SST) 4 x SDI (2 x SDI 3G), (2 x SDI 12G) 1 x DVI (single link, HDMI 1.4, HDCP 1.4)
Output Signal	1 x DP 1.2 (SST) 4 x SDI (2 x SDI 3G), (2 x SDI 12G) 1 x DVI (single link)
Power Supply	SMPS (AC 100 ~240V)
Power Consumption	350W max, 160W typ
Accessories Included	User manual, AC power cord, DP cable, HDMI cable, DVI-D cable, BNC cable, 3D glasses, Remote control and batteries
Unit Dimension	1268.5(W) x 753.3(H) x 85(D) mm 49.94(W) x 29.66(H) x 3.35(D) inch
Package Dimension	1453(W) x 937(H) x 300(D) mm 57.2(W) x 36.89(H) x 11.8(D) inch
Weight	33 kg, 72.75 lbs. (monitor) 43.8 kg, 96.56 lbs. (shipping package)
Certifications	ANSI/AAMI ES60601-1(UL60601-1), CAN/CSA C22.2 No. 60601-1, FCC Class B, EN60601-1, EN60601-1-2, CE, MDR(EU) 2017/745 Class I Medical Device

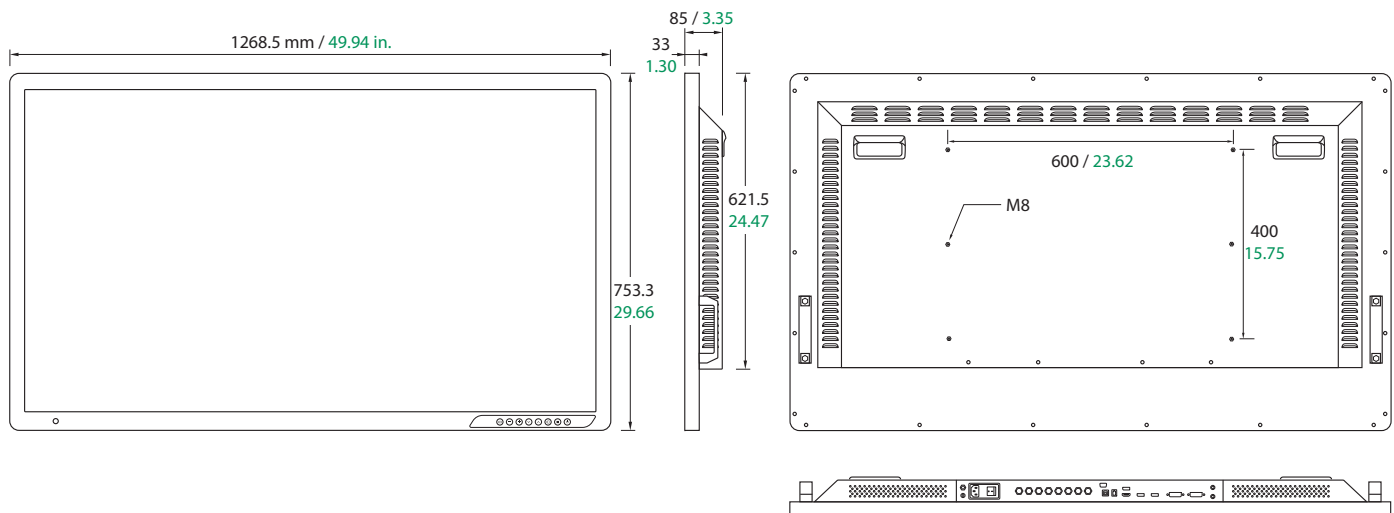

fsnmed.com

Specifications are subject to change with or without notice. Doc. # FSN2032 Rev. 3/25

Foreseeson Custom Displays, Inc.

2210 E. Winston Road
Anaheim, CA 92806 USA
Tel: 714-300-0540
Fax: 714-300-0546

Latin America
Aventura, FL 33180 USA
Tel: 714-507-3855

Foreseeson Korea

B-408, U-Space2
670 Daewangpangyo-ro,
Bundang-gu, Seongnam-si
Gyeonggi-do, Republic of Korea
Tel: +82-31-8017-0780

Foreseeson GmbH

Industriestrasse 38a
63150 Heusenstamm, Germany
Tel: +49 6104 64398 0
Fax: +49 6104 64398 11

Foreseeson UK Ltd.

1 Wolsey Road
East Molesay
Surrey
KT8 9EL
United Kingdom
Tel: +44 (0) 208 546 1047

Foreseeson (Shanghai) Medical Equipment Co., Ltd.

Room 1010, Building A
1439 Wuzhong Road
Rhein Hongjing Center
Minhang District, Shanghai, China
Tel: 18521095596